



## Rüsch CrystalClear Plus Tracheostomy Tube

### Supports Effective Tracheostomy Care

The single-use Rüsch CrystalClear Plus Tracheostomy Tube portfolio from Teleflex is comprised of cuffed, uncuffed, or fenestrated (cuffed and uncuffed) tracheostomy tubes, which are sold with two single-use, flexible inner cannulas. The range is ideally suited for a new Tracheostomy stoma, 'first tube in' situations. Manufactured from thermo-sensitive PVC tubing, and with a thin, flexible, PVC flange, Rüsch CrystalClear Plus Tracheostomy Tubes have been designed for patient comfort and effective stoma care<sup>1</sup>. The flexibility of the flange enables it to conform to the anatomy with minimal skin contact, when a stoma dressing is not required; it also helps to minimise tube movement during care, and provides ease of access to the stoma site for cleaning, to help reduce the risk of infection<sup>2</sup>.

### Designed with Patient Safety in Mind

Rüsch CrystalClear Plus Tracheostomy Tubes have size-specific colour-coded flexible inner cannulas that correspond to coloured writing on the flange. This helps to avoid the inadvertent use of an inappropriately sized inner cannula, which could increase 'work of breathing' and potentially compromise patient safety. Moreover, the flexible inner cannula can be removed and replaced, helping to prevent the accumulation of biofilm and secretions, which could occlude the outer cannula and restrict ventilation. Additional safety features include a swivel on the 15 mm connector to reduce torque on the ventilator circuit, and the 95° curvature of the tube, which facilitates a secure fit in new stomas to reduce the risk of accidental decannulation.

### Benefits



**The Clinician**  
Colour-coded inner cannulas aid size identification, for increased peace of mind



**Your Institution**  
Designed to reduce the risk of inadvertent tube displacement in the intensive care unit



**The Patient**  
Design features assist in maintaining a patent airway

1. Wilkinson KA, Martin IC, Freeth H, Kelly K, Mason M. NCEPOD: *On the right Trach?* 2014. Available at: <http://bit.ly/ontherighttrach>. Accessed August 22, 2017.  
2. Mallett J, Albarran JW, Richardson A, eds. *Critical Care Manual of Clinical Procedures and Competencies*. Oxford, UK: Wiley Blackwell; 2013.


# Rüsch CrystalClear Plus Tracheostomy Tube

## 15 mm Swivel Connector

Designed to reduce torque with ventilator circuits

## Cuff

Tubes are available both with or without a cuff

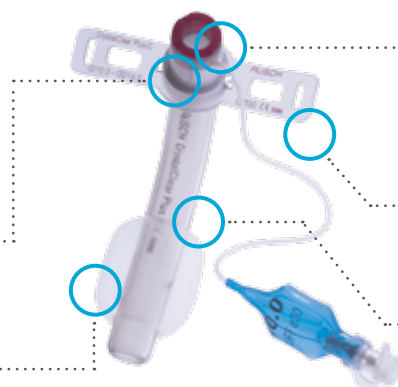

## Inner Cannula

Size-specific colour coding to aid size identification

## Thin, Flexible Flange

Facilitates effective stoma care

## 95° Shaft

Promotes a secure fit to help prevent accidental decannulation

## Rüsch CrystalClear Plus Tracheostomy Tube: Fast Facts

Single-use portfolio to help reduce the risk of infection associated with reusable tracheostomy tubes.

Sold with two flexible, colour-coded inner cannulas.

Fenestrated tracheostomy tubes available to support ventilator weaning and speech.

Compatible with the Rüsch PercuTwist Percutaneous Tracheostomy Dilator.

## Rüsch CrystalClear Plus Tracheostomy Tube

PRODUCT CODE*	PRODUCT CODE**	PRODUCT CODE*	PRODUCT CODE**	SIZE / I.D. (MM)	O.D. (MM)	LENGTH (MM)	ANGLE Ø APPROX.	CUFF Ø (MM)	INNER CANNULA	COLOUR CODE	QTY
UN-FENESTRATED				FENESTRATED							
121300-000035	121303-000035			3.5	5.3	47.0	95°	11.0	-	-	5
121300-000040	121303-000040			4.0	5.9	54.0	95°	12.0	-	-	5
121300-000045	121303-000045			4.5	6.5	54.0	95°	14.0	-	-	5
121300-000050	121303-000050			5.0	7.0	57.0	95°	14.0	-	-	5
121300-000055	121303-000055			5.5	7.5	58.0	95°	15.0	-	-	5
121300-000060	121303-000060	121301-000060	121302-000060	6.0	8.7	64.0	95°	19.0	6.0	Orange	5
121300-000065	121303-000065	121301-000065	121302-000065	6.5	8.7	64.0	95°	19.0	6.5	Maroon	5
121300-000070	121303-000070	121301-000070	121302-000070	7.0	10.0	70.0	95°	21.0	7.0	Blue	5
121300-000075	121303-000075	121301-000075	121302-000075	7.5	10.7	75.5	95°	25.0	7.5	Yellow	5
121300-000080	121303-000080	121301-000080	121302-000080	8.0	11.3	81.0	95°	28.0	8.0	Lilac	5
121300-000085	121303-000085	121301-000085	121302-000085	8.5	11.3	81.0	95°	28.0	8.5	Green	5
121300-000090	121303-000090	121301-000090	121302-000090	9.0	12.7	87.0	95°	28.0	9.0	Brown	5
121300-000095	121303-000095	121301-000095	121302-000095	9.5	12.7	87.0	95°	28.0	9.0	Navy	5
121300-000100	121303-000100	121301-000100	121302-000100	10.0	14.0	95.5	95°	30.0	10.0	Red	5

\* Without Cuff. \*\* With Cuff.

## Rüsch CrystalClear PDT Set

PRODUCT CODE	SIZE / I.D. (MM)	O.D. (MM)	LENGTH (MM)	ANGLE Ø APPROX.	CUFF Ø (MM)	COLOUR CODE	QTY
UN-FENESTRATED WITH CUFF							
121304-000070	7.0	10.0	70.0	95°	21.0	Blue	1
121304-000080	8.0	11.3	81.0	95°	28.0	Lilac	1
121304-000090	9.0	12.7	87.0	95°	28.0	Brown	1

## Rüsch CrystalClear Accessories

PRODUCT CODE	DESCRIPTION	QTY
121311	Shower cap for Rüsch CrystalClear Plus	4
121312	Cough cap for Rüsch CrystalClear Plus	4
121313	Sealing cap for Rüsch CrystalClear Plus	4

## Rüsch CrystalClear Inner Cannula

PRODUCT CODE	SIZE	I.D. (MM)	O.D. (MM)	COLOUR CODE	QTY
121310-000060	6.0	5.0	5.7	Orange	10
121310-000065	6.5	5.0	5.7	Maroon	10
121310-000070	7.0	6.0	6.7	Blue	10
121310-000075	7.5	6.5	7.3	Yellow	10
121310-000080	8.0	7.0	7.7	Lilac	10
121310-000085	8.5	7.0	7.7	Green	10
121310-000090	9.0	8.0	8.7	Brown	10
121310-000095	9.5	8.0	8.7	Navy	10
121310-000010	10.0	8.8	9.7	Red	10



Contents of unopened, undamaged package are clean, ready to use and for single-use. Refer to package insert for current warnings, indications, contraindications, precautions, and Instructions for Use.

Teleflex, the Teleflex logo, CrystalClear, PercuTwist, and Rüsch are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative.

Revised: 08/2017 © 2017 Teleflex Incorporated. All rights reserved. 94 30 67 - 00 00 01 - REV A - 08 17 00.25



Distributed by:  
Teleflex Headquarters International, Ireland · Teleflex Medical Europe Ltd. · IDA Business & Technology Park  
Dublin Road · Athlone · Co Westmeath · Tel. +353 (0)9 06 46 08 00 · Fax +353 (0)14 37 07 73 · orders.intl@teleflex.com  
United Kingdom Tel. +44 (0)14 94 53 27 61 · info.uk@teleflex.com  
South Africa Tel. +27 (0)11 807 4887 · assist.africa@teleflex.com