



**Hudson RCI**  
**Softech Plus / ETCO<sub>2</sub> Sampling Cannula**  
 For reliable oxygen delivery and ETCO<sub>2</sub> sampling

**Confidence in oxygen delivery**


At Teleflex, we continuously strive to improve the standard of safety and care for patients. Our innovative Hudson RCI Nasal Cannula Solutions are designed to meet the challenges presented by unique patient needs and physician preferences.

**Softech Plus ETCO<sub>2</sub> Sampling Cannula**

The Softech Plus ETCO<sub>2</sub> Sampling Cannula from Teleflex uses a divided design to deliver oxygen from one nare and sample carbon dioxide from the other, permitting an undiluted end tidal gas sample.

Clinicians should feel confident in choosing Softech Plus Cannulas for comfortable oxygen delivery because of the proprietary combination of a soft, non-DEHP material and improved product design. Softech Plus Cannulas can help improve patient comfort and enhance quality of care.

The Softech Plus Material is incorporated into the nasal prongs and entire bolo-style lariat, combining comfort with an ideal fit for most patients.

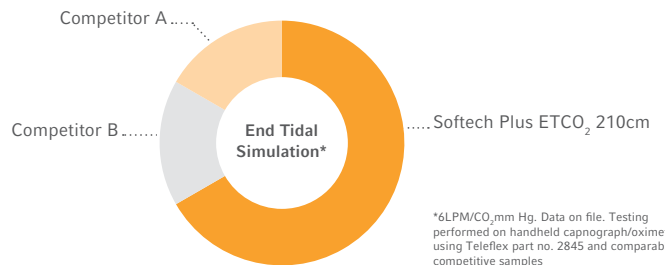

**Benefits**




**The Clinician**  
 Split nare design helps produce a sharp wave form for clinician diagnostics



**Your Institution**  
 Range of connection options to fit clinical/facility needs



**The Patient**  
 Softer nasal prongs and lariat

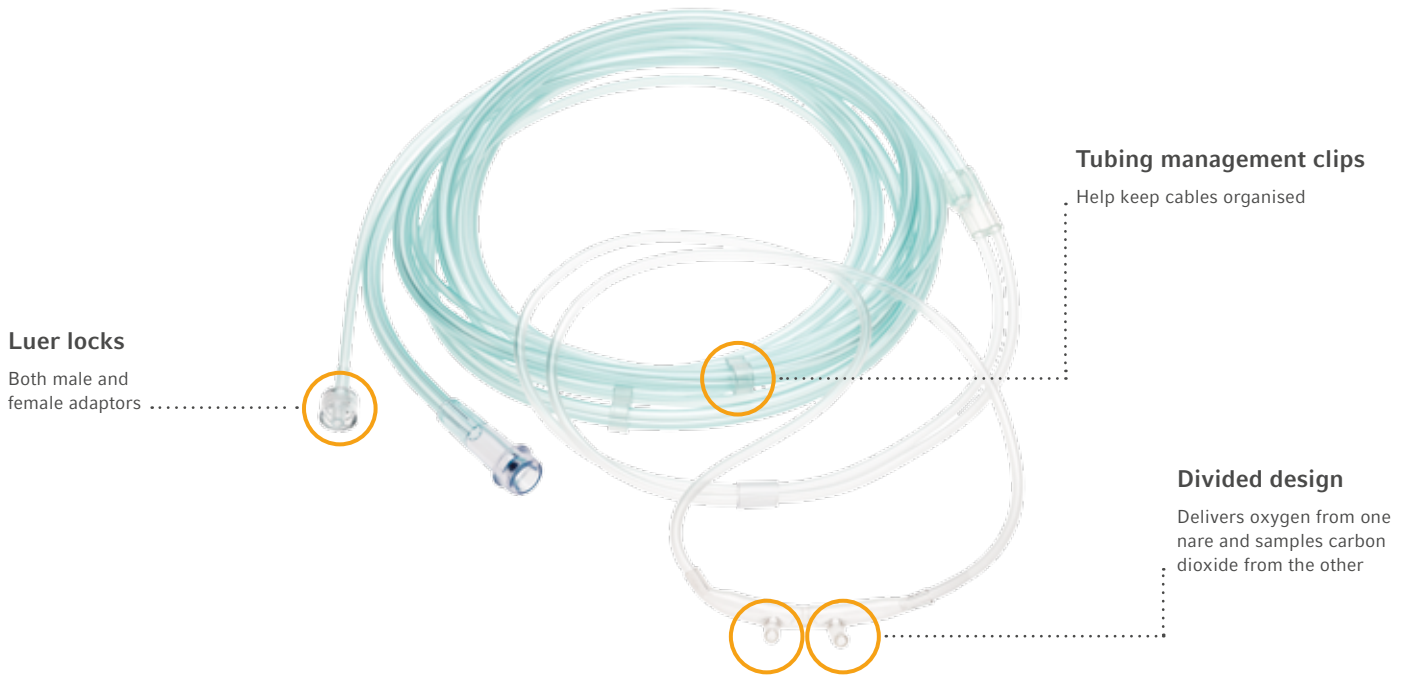

**Luer locks**

Both male and female adaptors

**Tubing management clips**

Help keep cables organised

**Divided design**

Delivers oxygen from one nare and samples carbon dioxide from the other

**Hudson RCI Softech Plus ETCO<sub>2</sub> Sampling Cannula Fast Facts**

The Softech Plus ETCO<sub>2</sub> Sampling Cannula features softer nasal prongs and lariat for greater patient comfort

The non-DEHP construction of Softech Plus cannulas is not made with natural rubber latex

Smaller nare design allows for enhanced patient comfort. Now available in adult, paediatric and infant sizes

Helps enhance clinical efficiency with a range of connection options to fit facility needs, while tubing management clips help keep oxygen delivery organised

Helps enhance clinical efficiency with light green tubing and pigtails that vary from 50 mm to 420 cm to connect to different types of machines

**Hudson RCI Softech Plus ETCO<sub>2</sub> Sampling Cannula**

PART NO.	CANNULA DESCRIPTION	LENGTH OF O <sub>2</sub> SAMPLING LINE	LENGTH OF CO <sub>2</sub> SAMPLING LINE - PIGTAIL	LUER CONNECTOR
2842	Adult Divided	7' (210cm)	2" (50.8mm)	Male
2843	Adult Divided	7' (210cm)	2" (50.8mm)	Female
2844	Adult Divided	7' (210cm)	7' (210cm)	Female
2845	Adult Divided	7' (210cm)	7' (210cm)	Male
2846	Adult Divided	14' (420 cm)	2" (50.8mm)	Female
2847	Adult Divided	14' (420 cm)	14' (420 cm)	Female
2848	Adult Divided	14' (420 cm)	14' (420 cm)	Male
2849	Paediatric Divided	7' (210cm)	2" (50.8mm)	Female
2850	Paediatric Divided	7' (210cm)	7' (210cm)	Female
2851	Paediatric Divided	7' (210cm)	7' (210cm)	Male
2852	Paediatric Divided	14' (420 cm)	2" (50.8mm)	Female
2853	Paediatric Divided	14' (420 cm)	14' (420 cm)	Female
2854	Paediatric Divided	14' (420 cm)	14' (420 cm)	Male
2855	Infant Divided	7' (210cm)	2" (50.8mm)	Female
2856	Infant Divided	7' (210cm)	7' (210cm)	Female
2857	Infant Divided	7' (210cm)	7' (210cm)	Male



**Softech Plus ETCO<sub>2</sub> Nares**

The Softech Plus ETCO<sub>2</sub> Sampling Cannula delivers oxygen from one Nare and samples expired gases from the other.

940827 - 00 00 01 · REV A · 01 / 16 · 02 16 xx Teleflex, the Teleflex logo, Hudson RCI and Softech are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Information in this document is not a substitute for the product Instructions for Use. The products in this catalogue may not be available in all countries. Please contact your local representative. All data current at time of printing (01/2016). Subject to technical changes without further notice. © 2016 Teleflex Incorporated. All rights reserved.



Distributed by:  
**Teleflex Headquarters International**, Ireland · Teleflex Medical Europe Ltd. · IDA Business & Technology Park  
 Dublin Road · Athlone · Co Westmeath · Tel. +353 (0)9 06 46 08 00 · Fax +353 (0)14 37 07 73 · orders.intl@teleflex.com  
 United Kingdom Tel. +44 (0)14 94 53 27 61 · info.uk@teleflex.com