



ARROW PAEDIATRIC CVC SETS

Nitinol guidewire and safety features



TELEFLEX'S ARROW PAEDIATRIC CVC SETS: THE FIRST CHOICE FOR VENOUS ACCESS CVC CATHETER NOW COMES WITH ADVANCED COMPONENTS


Although generally safe and routine in ICU, central venous catheter insertion procedures can be challenging and involve complications, especially where a child is concerned.

Teleflex's ARROW has always strived to provide the optimum quality products to ensure ease-of-use and the utmost in patient safety and comfort. That is why ARROW central venous catheters are the first choice for venous access. ARROW's paediatric CVCs are, of course, no exception. The flexi-

ble polyurethane they are made of is durable but pliant and soft enough to protect the vessel.



ARROW's paediatric CVCs come in multiple lumen and length configurations to meet the specific therapy needs of children. Available configurations range from 1 to 3 lumens and come in several lengths and French sizes. You also have the choice of uncoated or coated catheters for antimicrobial protection.

TELEFLEX'S ARROW PAEDIATRIC CVC SETS WITH NITINOL GUIDEWIRE, INTRODUCER NEEDLE AND ADDITIONAL SAFETY FEATURES

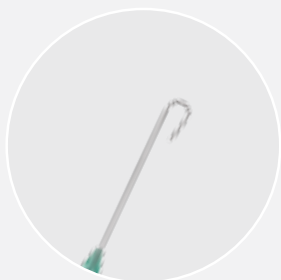

In order to meet customer needs, we have reviewed the content of our paediatric CVC sets and now offer the following:

- Nitinol guidewire: the Nitinol guidewire has a welded connection at the tip, a high tensile strength and is highly flexible. Additionally the guidewire has cm-markings along the shaft and a small J tip radius.
- New introducer needle: the introducer needle provided is siliconised and has a sharper bevel, this introducer is used in our other products. These features can provide more comfort for your patient and make it easier for the clinical staff to attend to the task at hand.
- SharpsAway II Locking Disposal Cup and safety scalpel: the avoidance of needle stick injuries for your staff not only complies with key regulatory guidelines, it's being responsible by minimising the risk to your patient and the clinical staff delivering life-saving care. Our safety scalpel is easy-to-control to avoid sharps injuries, while our SharpsAway II® Locking Disposal Cup allows safe disposal of used sharps to minimise the risk of exposure.

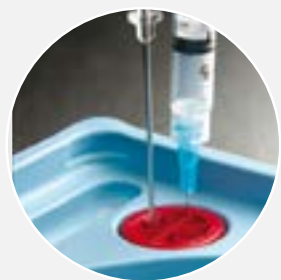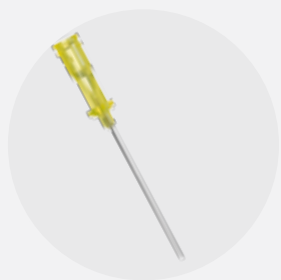

ANTIMICROBIAL TECHNOLOGY
(SEE LIST OF CODES BELOW
FOR AGB VERSIONS)

patented treatment bonds chlorhexidine and silver sulfadiazine to the entire surface length



BLUE FLEXTIP®

ARROW's soft, contoured tip is designed to be easier to manoeuvre and minimise trauma

MOULDED, FLEXIBLE, LOW-PROFILE JUNCTION HUB

designed for patient comfort and ease of maintenance

SET COMPONENTS

- polyurethane catheter (PUR) with Blue FlexTip
- SWG Nitinol guidewire and ARROW Advancer™
- introducer catheter over needle
- 5 ml Luer slip syringe
- introducer needles
- 2 dilators
- catheter clamp and fastener
- safety scalpel
- SharpsAway II Locking Disposal Cup



FLEXIBLE POLYURETHANE

designed to soften in situ and enhance indwelling characteristics

PAEDIATRIC NITINOL GUIDEWIRE KITS – WITH CM-MARKING, LATEX-FREE



REF.	CATHETER	CROSS SECTION	LUMEN	DIAMETER	TIP	MATERIAL	QTY
EU-04522-N	22 Ga x 4" (10 cm)	n/a	1	.021" (0.53 mm) x 13¾" (35 cm)	Straight/Straight Nitinol	PUR	5
EU-04150-N	20 Ga x 5" (12 cm)	n/a	1	.025" (0.64 mm) x 13¾" (35 cm)	J/straight Nitinol	PUR	5
EU-04218-N	18 Ga x 8" (20 cm)	n/a	1	.018" (0.46 mm) x 17¾" (45 cm)	J/straight Nitinol	PUR	5
EU-12402-N	4 Fr x 2" (5 cm)	C-22/22	2	.018" (0.46 mm) x 17¾" (45 cm)	J/straight Nitinol	PUR	5
EU-15402-N	4 Fr x 3⅞" (8 cm)	C-22/22	2	.018" (0.46 mm) x 17¾" (45 cm)	J/straight Nitinol	PUR	5
EU-25402-N	4 Fr x 3⅞" (8 cm)	C-22/22	2	.018" (0.46 mm) x 17¾" (45 cm)	J/straight Nitinol	PUR/AGB	5
EU-14402-N	4 Fr x 5⅞" (13 cm)	C-22/22	2	.018" (0.46 mm) x 17¾" (45 cm)	J/straight Nitinol	PUR	5
EU-24402-N	4 Fr x 5⅞" (13 cm)	C-22/22	2	.018" (0.46 mm) x 17¾" (45 cm)	J/straight Nitinol	PUR/AGB	5
EU-16402-N	4 Fr x 12" (30 cm)	C-22/22	2	.018" (0.46 mm) x 26¾" (68 cm)	J/straight Nitinol	PUR	5
EU-14502-N	5 Fr x 5⅞" (13 cm)	E-18/20	2	.021" (0.53 mm) x 17¾" (45 cm)	J/straight Nitinol	PUR	5
EU-22502-N	5 Fr x 2" (5 cm)	E-18/20	2	.021" (0.53 mm) x 17¾" (45 cm)	J/straight Nitinol	PUR/AGB	5
EU-15403-N	4 Fr x 3" (8 cm)	H-20/22/22	3	.018" (0.46 mm) x 17¾" (45 cm)	J/straight Nitinol	PUR	5
EU-14403-N	4 Fr x 5⅞" (13 cm)	H-20/23/23	3	.018" (0.46 mm) x 17¾" (45 cm)	J/straight Nitinol	PUR	5
EU-15553-N	5,5 Fr x 3" (8 cm)	H-20/22/22	3	.018" (0.46 mm) x 17¾" (45 cm)	J/straight Nitinol	PUR	5
EU-25553-N	5,5 Fr x 3⅞" (8 cm)	H-20/22/22	3	.018" (0.46 mm) x 17¾" (45 cm)	J/straight Nitinol	PUR/AGB	5
EU-16553-N	5,5 Fr x 5⅞" (13 cm)	H-20/22/22	3	.018" (0.46 mm) x 17¾" (45 cm)	J/straight Nitinol	PUR	5
EU-26553-N	5,5 Fr x 5⅞" (13 cm)	H-20/22/22	3	.018" (0.46 mm) x 17¾" (45 cm)	J/straight Nitinol	PUR/AGB	5
EU-14553-N	5,5 Fr x 12" (30 cm)	H-20/22/22	3	.018" (0.46 mm) x 26¾" (68 cm)	J/straight Nitinol	PUR	5

AGB = ARROWg*ard Blue

Teleflex is a leading global provider of specialty medical devices used for diagnostic and therapeutic procedures in critical care, urology and surgery. Our mission is to provide solutions that enable healthcare providers to improve outcomes and enhance patient and provider safety. We specialise in devices for general and regional anaesthesia, cardiac care, respiratory care, urology, vascular access and surgery and we serve healthcare providers in more than 150 countries. Teleflex also provides specialty products for medical device manufacturers.

Our well known brands include **ARROW**[®], **DEKNATEL**[®], **GIBECK**[®], **HUDSON RCI**[®], **KMEDIC**[®], **LMA**[™], **PILLING**[®], **PLEUR-EVAC**[®], **RÜSCH**[®], **SHERIDAN**[®], **TAUT**[®], **TFX OEM**[®], **VASONOVA**[™], **VIDACARE**[®] and **WECK**[®], all of which are trademarks or registered trademarks of Teleflex Incorporated.

Teleflex global operations: Australia, Austria, Belgium, Canada, China, Czech Republic, France, Germany, Greece, India, Ireland, Italy, Japan, Malaysia, Mexico, Netherlands, New Zealand, Portugal, Singapore, Slovak Republic, South Africa, Spain, Switzerland, United Kingdom, Uruguay and USA.

YOUR INTERNATIONAL CONTACTS:

TELEFLEX HEADQUARTERS INTERNATIONAL, IRELAND

Teleflex Medical Europe Ltd., IDA Business and Technology Park,
Dublin Road, Athlone, Co Westmeath
Phone +353 (0)9 06 46 08 00 · Fax +353 (0)14 37 07 73
orders.intl@teleflex.com

AUSTRALIA/NEW ZEALAND +61 (0)3 9081 0600

AUSTRIA +43 (0)1 402 47 72

BELGIUM +32 (0)2 333 24 60

CHINA (SHANGHAI) +86 (0)21 6163 0965

CHINA (BEIJING) +86 (0)10 6418 5699

CZECH REPUBLIC +420 (0)495 759 111

FRANCE +33 (0)5 62 18 79 40

GERMANY +49 (0)7151 406 0

GREECE +30 210 67 77 717

INDIA +91 (0)44-2836 5040

ITALY +39 0362 58 911

JAPAN +81 (0)3 3379 1511

NETHERLANDS +31 (0)88 00 215 00

PORTUGAL +351 22 541 90 85

SINGAPORE +65 6439 3000

SLOVAK REPUBLIC +421 (0)3377 254 28

SOUTH AFRICA +27 (0)11 807 4887

SPAIN +34 918 300 451

SWITZERLAND +41 (0)31 818 40 90

UNITED KINGDOM +44 (0)1494 53 27 61

For detailed information see www.teleflex.com

Arrow, SharpsAway and Advancer are trademarks or registered (in the U.S. and/or other countries) trademarks of Teleflex Incorporated or its affiliates.

The products in this catalogue may not be available in all countries. Please contact your local representative. All data current at time of printing (10/2014). Subject to technical changes without further notice.

94 11 26 - 00 00 01 · REV A · MC / WM · 10 14 01