

CATHETER MOUNTS



The ideal variety of products for connection of the breathing system or HME/filter to the patient

The range consists of three high quality flextubes, connected to different adaptors and a double swivel.


PRODUCT INFORMATION



Smooth-Flo® with double swivel

Swivel with sealaround 9 mm cap for fiber optic bronchoscopy and 3.5 mm port for suctioning

Connectors: patient side 22 M / 15 F
machine side 22 M

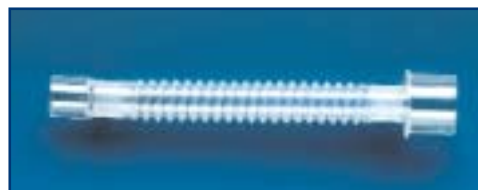

Smooth-Flo® Flex

packed with red dust cap for protection

Connectors: patient side 15 F
machine side 22 F

AT A GLANCE

Smooth-Flo®

Smoothbore type flex-tube minimizes water collection and therefore the growth of microorganisms.

Soft-Flex

Soft flexible tubing decreases the risk of constriction of the inner lumen through kinking or traction.

Expandi-Flex

Continuously expandable flextube that offers a variety of lengths and can be fixed in different positions (Memory Function).

PRODUCT INFORMATION



Soft-Flex with double swivel

Swivel with sealaround 9 mm cap for fiber optic bronchoscopy and 3.5 mm port for suctioning

Connectors: patient side 22 M / 15 F
machine side 22 F



Expandi-Flex with double swivel

Swivel with sealaround 9 mm cap for fiber optic bronchoscopy and 3.5 mm port for suctioning

Connectors: patient side 22 M / 15 F
machine side 22 M




Soft-Flex with elbow

Connectors: patient side 22 M / 15 F
machine side 22 F



Expandi-Flex with elbow

Connectors: patient side 22 M / 15 F
machine side 22 F

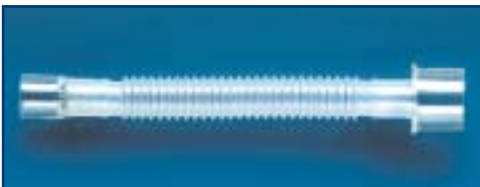

Soft-Flex

packed with red dust cap for protection,

Connectors: patient side 15 F
machine side 22 F



Expandi-Flex

packed with red dust cap for protection,

Connectors: patient side 15 F
machine side 22 F

ORDERING INFORMATION

Ref.	Description	Weight	Length	Dead space	Qty
22511	Smooth-Flo® with double swivel, sterile	26.1 g	205 mm	34 ml	25
22512	Smooth-Flo® with double swivel, clean	26.1 g	205 mm	34 ml	35
22521	Soft-Flex with double swivel, sterile	20.5 g	210 mm	40 ml	25
22522	Soft-Flex with double swivel, clean	20.5 g	210 mm	40 ml	35
22531	Expandi-Flex with double swivel, sterile	20.1 g	140/200 mm	23/39 ml	30
22532	Expandi-Flex with double swivel, clean	20.1 g	140/200 mm	23/39 ml	40
22542	Expandi-Flex with elbow, clean	16.8 g	125/180 mm	27/43 ml	40
22552	Smooth-Flo® Flex, clean	16.4 g	187 mm	38.5 ml	35
22562	Expandi-Flex, clean	14.3 g	122/179 mm	22/38 ml	25
22572	Soft-Flex with elbow, clean	17.3 g	195 mm	42.5 ml	35
22582	Soft-Flex, clean	14.7 g	195 mm	36 ml	35
23411	Double swivel connector, sterile	9.7 g	45 mm	6.5 ml	20
23412	Double swivel connector, clean	9.7 g	45 mm	6.5 ml	20

AT A GLANCE

- high quality see-through material
- 22 female connector machine side fits over 22 mm male Y-piece
- can be combined with Gibeck HMEs & bacterial/viral filters
- environmental-friendly after incineration
- latex-free

ACCESSORIES



Double swivel connector

Swivel with sealaround 9 mm cap for fiber optic bronchoscopy and 3.5 mm port for suctioning,

Connectors: patient side 22 M / 15 F
machine side 22 M