

Arrow  
**Epidural Kit and CSE Kit**  
 The Quality of Care You Demand for Your Patients

### Control Complications in the Epidural Space

Teleflex continuously seeks to improve techniques and technologies in post-operative and long-term pain management. Regional anesthesia, which allows for controlled patient stress, is gaining importance. It helps to lessen the risk of heart and lung complications while reducing operating theatres waiting times.

### Reliable Performance

The Arrow Epidural Kit from Teleflex is composed of a Tuohy Epidural Needle meticulously constructed to make it reliable, and a nylon epidural catheter with a multiport tip configuration. Anaesthetic is injected through the catheter and is dispersed through the three ports located at its distal end. The catheter is designed to be kink-resistant and to have a rounded and smooth

tip to control the risk of trauma. Additionally, the catheter is designed to be rigid enough to allow insertion through the needle into the epidural space without the use of a stylet. A radiopaque stripe ensures visibility when needed.

### Procedure in Combined Epidural/Spinal Anaesthesia

The Arrow CSE (combined spinal-epidural) Kit from Teleflex has been developed to help you perform a CSE procedure and meet the needs of your patients. It allows patients to experience less trauma by benefiting from a single puncture when needed. The rapid onset of spinal anesthesia reduces operating theatres waiting time.

#### Benefits




**The Clinician**  
 Shorter procedural time.



**Your Institution**  
 Effective pain management reduces length of stay in hospital.



**The Patient**  
 Less Trauma.  
 Effective pain management.  
 Faster recovery.

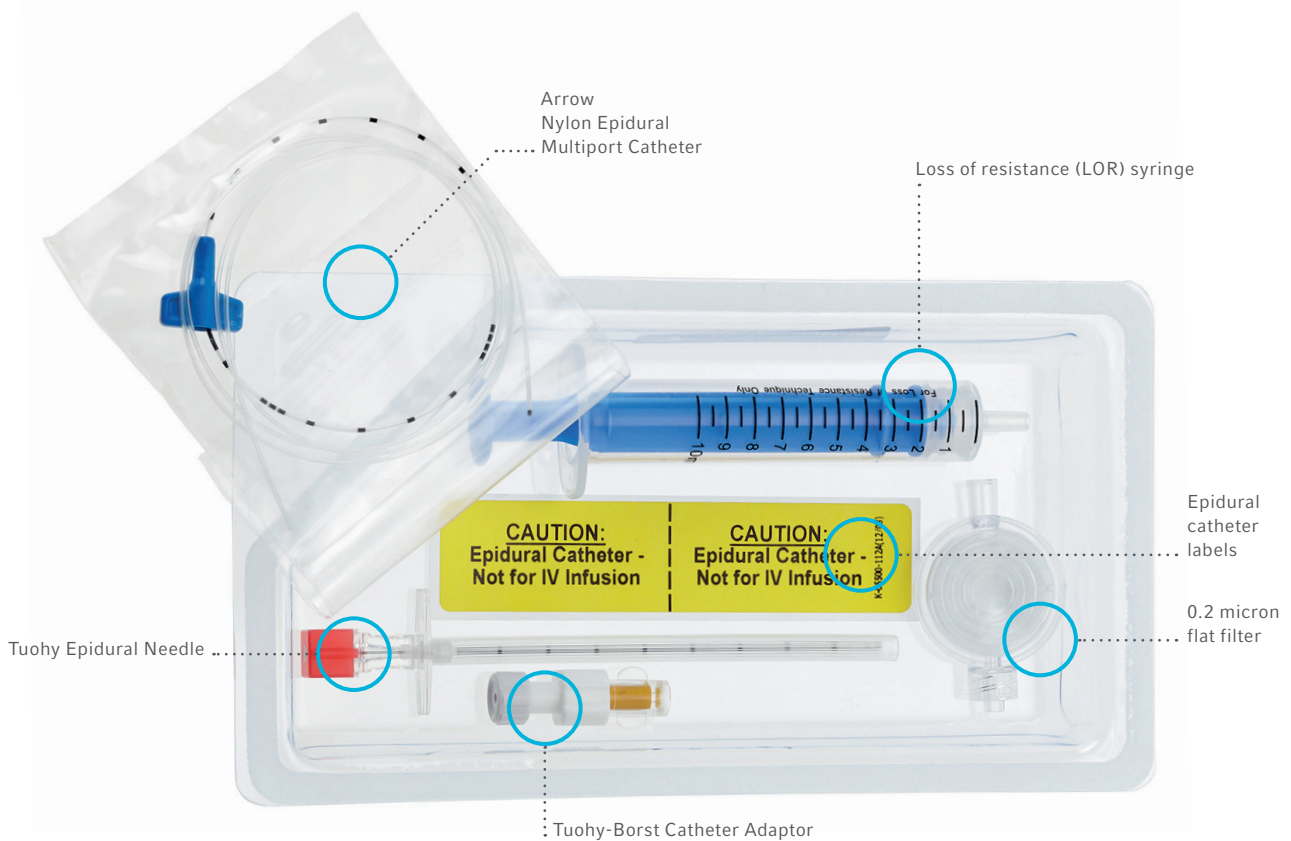

## Epidural Kit and CSE Kit Fast Facts

Thread assist device to facilitate catheter insertion by allowing for precision and versatility

Epidural catheter with 3 distal ports evenly positioned around the catheter body ensuring easy fluid aspiration

Pencil point spinal needle in CSE kit enables smooth insertion

Sterile  
Not made with natural rubber latex  
Not made with DEHP

### Epidural Kit and Arrow CSE Kit:

REF.	LOR SYRINGE 10 ML	SPINAL NEEDLE PENCIL POINT 27 G X 124 MM	UNITS PER CASE
<b>Epidural Set</b>			
MPN-18020-S	x		20
MPN-18020			20
<b>CSE Set</b>			
MPC-18020-S	x	x	20
MPC-18020		x	20

### Set Components:

- 1 x 20 G x 90 cm Arrow Nylon Epidural Multiport Catheter with centimeter markings and thread assist device
- 1 x 18 G x 90 mm Tuohy Epidural Needle with fixed wings and centimeter markings
- 1 x 27 G Pencil Point Spinal Needle (MPC-18020-S & MPC-18020 only)
- 1 x 0.2 micron flat filter
- 1 x Tuohy-Borst Catheter Adaptor
- 2 x Epidural catheter labels

Teleflex, the Teleflex logo and Arrow are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Information in this document is not a substitute for the product Instructions for Use. The products in this document may not be available in all countries. Please contact your local representative. All data current at time of printing (02/2016). Subject to technical changes without further notice. © 2016 Teleflex Incorporated. All rights reserved. 94 08 00 - 00 00 01 - REV A - 02 16 01