



More Radial,
Less Restrictions

GuideLiner[®] V3 Catheter
5F Guide Extension Catheter

End-to-End Innovation Supports the Radial Reach

True Support for Small Vessel Delivery

The 5F GuideLiner® V3 Catheter was designed specifically to support 5F guide interventions and small vessel delivery. Its small crossing profile allows it to fit through any 5F guide catheter, which means fewer guide catheter restrictions and greater confidence during interventions.



Increases back-up support while allowing deep-seating for coaxial alignment, facilitation of distal device delivery and selective delivery of contrast



0.046" I.D. to Facilitate Select Device Delivery*

3mm Stent Crossing Profiles

MANUFACTURER	STENT MODEL	DIAMETER	AVERAGE MEASURED CROSSING PROFILE**	MAXIMUM MEASURED CROSSING PROFILE**
Boston Scientific	Ion™	3 mm	0.042"	0.044"
Boston Scientific	PROMUS Element™ Plus	3 mm	0.041"	0.042"
Boston Scientific	VeriFLEX™	3 mm	0.040"	0.040"
Abbott	Xience Prime® LL	3 mm	0.043"	0.043"


*Stent diameters >3mm are not recommended for delivery in conjunction with the 5F GuideLiner® V3 Catheter.
The Twin-Pass® Dual Access Catheters are not compatible with the 5F GuideLiner® V3 Catheter.

**Data on file.


GuideLiner® catheters are intended to be used in conjunction with guide catheters to access discrete regions of the coronary and/or peripheral vasculature, and to facilitate placement of interventional devices.

MODEL	SIZE	REQUIRED GUIDE CATHETER I.D.	GIDELINER I.D.	RAPID EXCHANGE LENGTH	WORKING LENGTH
5569	5F	5F I.D. $\geq 0.056"$ (1.42mm)	0.046" (1.17mm)	25 cm	150cm

Packaged in quantities of 1 unit per box.

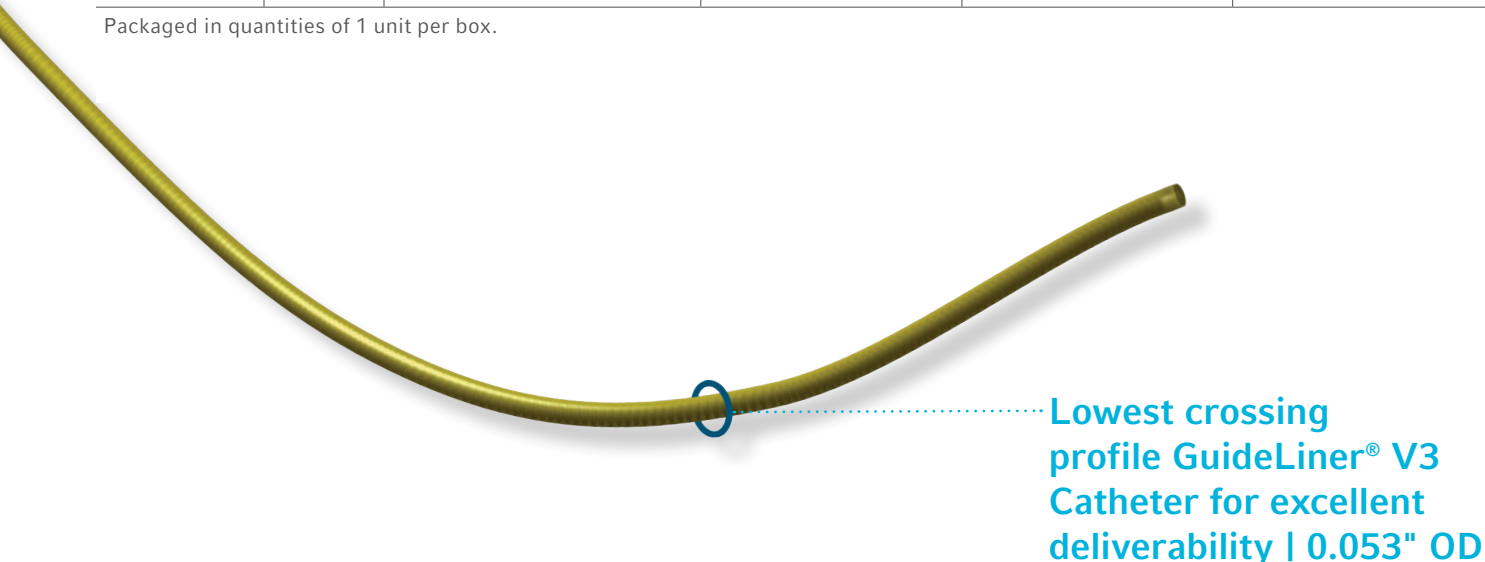

5F GuideLiner® V3 Catheter Dimensions

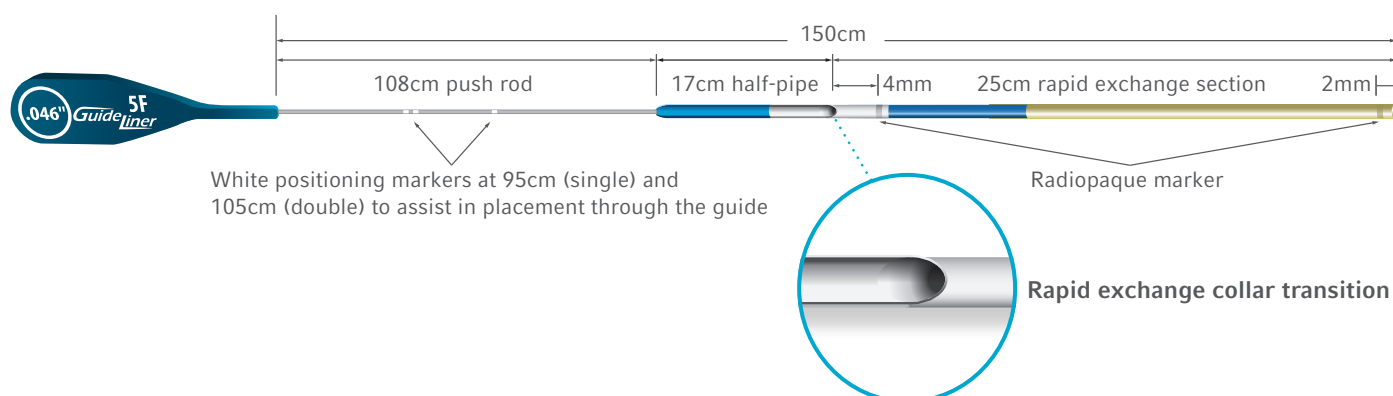

Teleflex is a global provider of medical technologies designed to improve the health and quality of people's lives. We apply purpose-driven innovation – a relentless pursuit of identifying unmet clinical needs – to benefit patients and healthcare providers. Our portfolio is diverse, with solutions in the fields of vascular and interventional access, surgical, anesthesia, cardiac care, urology, emergency medicine and respiratory care. Teleflex employees worldwide are united in the understanding that what we do every day makes a difference. For more information, please visit teleflex.com.

Teleflex is the home of Arrow®, Deknatel®, Hudson RCI®, LMA®, Pilling®, Rüsch®, UroLift®, and Weck® – trusted brands united by a common sense of purpose.

Corporate Office

Phone +1 610 225 6800, 550 E. Swedesford Road, Suite 400, Wayne, PA 19087, USA

Regional Offices

United States: Phone +1 919 544 8000, Toll Free 866 246 6990, cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

Latin America: Phone +1 919 433 4999, la.cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

International: Phone +353 (0)9 06 46 08 00, orders.intl@teleflex.com, Teleflex Medical Europe Ltd., IDA Business and Technology Park, Dublin Road, Athlone, Co Westmeath, Ireland

Australia 1300 360 226

Austria +43 (0)1 402 47 72

Belgium +32 (0)2 333 24 60

Canada +1 (0)800 387 9699

China (Shanghai) +86 (0)21 6163 0965

China (Beijing) +86 (0)10 6418 5699

Czech Republic +420 (0)495 759 111

France +33 (0)5 62 18 79 40

Germany +49 (0)711 2090 8000

Greece +30 210 67 77 717

India +91 (0)44 2836 5040

Italy +39 0362 58911

Japan +81 (0)3 6632 3600

Korea +82 2 536 7550

Mexico +52 55 5002 3500

Netherlands +31 (0)88 00 215 00

New Zealand 0800 601 100

Poland +48 22 4624032

Portugal +351 22 541 90 85

Singapore (SEA non-direct sales countries) +65 6439 3000

Slovak Republic +421 (0)3377 254 28

South Africa +27 (0)11 807 4887

Spain +34 918 300 451

Switzerland +41 (0)31 818 40 90

United Kingdom +44 (0)1494 53 27 61

For more information, please visit teleflex.com.

Please see the Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions

CAUTION: Federal (USA) law restricts this device to sale by or on the order of a physician.

Teleflex, the Teleflex logo, GuideLiner and Twin-Pass are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Other trademarks are the property of their respective owners.

Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative.

© 2019 Teleflex Incorporated. All rights reserved. MC-005671 LA EN