


PRONTO[®] LP

extraction catheter

The Pronto LP extraction catheter is indicated for the removal of fresh, soft emboli and thrombi from vessels in the coronary and peripheral system

Description	Model	Outer Diameter	Crossing Profile	Proximal Shaft Shape	Guide Catheter Compatibility	Extraction Lumen Area	Minimum Vessel Diameter	Box Quantity
Pronto LP	5010	0.059" (1.50mm)	0.053"	Round	6F (≥0.066")	0.69mm ²	1.50mm	1
Pronto LP 5F	5015	0.054" (1.37mm)	0.053"	U-Channel	5F (≥0.056")	0.69mm ²	1.50mm	1

Packaged in quantities of 1 unit per box.
The Pronto LP comes complete with two locking 30ml aspiration syringes and a 70µm filter basket.

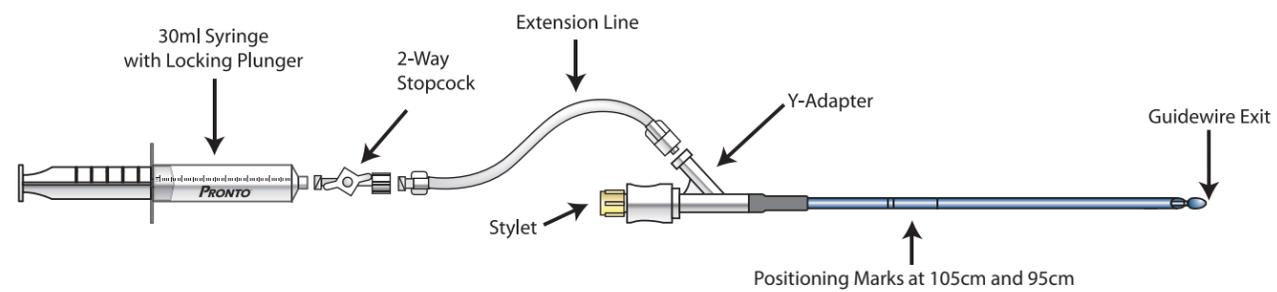

PRONTO[®] LP

extraction catheter

LOW-PROFILE ASPIRATION
now with a 5F guide compatibility option

CAUTION: Federal Law (U.S.A.) restricts this device to sale by or on the order of a physician.

Please see the Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions.

Pronto is a registered trademark of Vascular Solutions, Inc.
©2015 Vascular Solutions, Inc. All rights reserved. ML1873 Rev. E 12/15



Vascular Solutions, Inc.
6464 Sycamore Court North
Minneapolis, Minnesota 55369 USA
USA Customer Service: 888.240.6001
www.vasc.com



ADVANCED DESIGN IN ASPIRATION

The 0.014" compatible Pronto LP combines the extraction power of the Pronto in a low, 0.053" crossing profile for use in vessels $\geq 1.5\text{mm}$ in diameter. Two versions are available offering $\geq 6\text{F}$ or $\geq 5\text{F}$ guide compatibility for different clinical applications.

PATENTED SILVA TIP

provides a tapered entry profile that combines vessel protection with the ability to track through tortuous anatomy



SLOPED EXTRACTION LUMEN

accommodates larger thrombi


DISTAL HYDROPHILIC COATING

for lubricious delivery and smooth navigation



PRELOADED STYLET


can easily be removed for aspiration through a separate port, and offers pushability and kink resistance during delivery

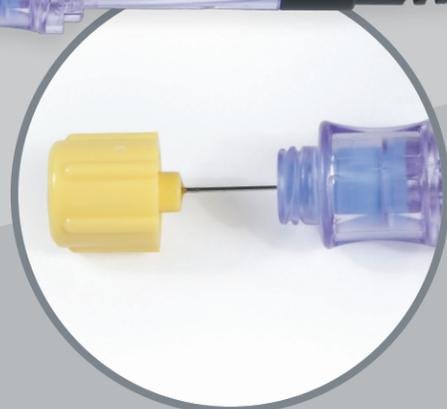


LOW-PROFILE OPTIONS WITH MAXIMUM PERFORMANCE

Original Pronto LP (Model 5010)

- Allows for the delivery of a second guidewire or contrast around the LP in a $\geq 6\text{F}$ guide catheter
- For delivery into distal, tortuous vessels
- Round, braided shaft and stylet for kink resistance

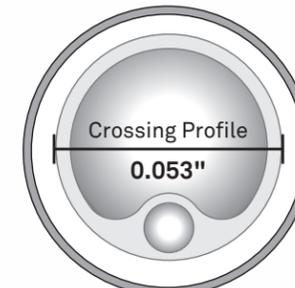

DISTAL SHAFT

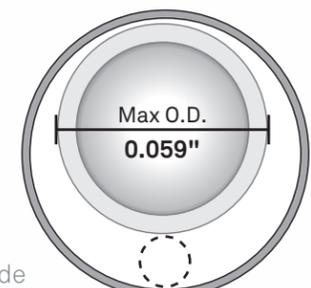

PROXIMAL SHAFT

shown in 6F guide

0.014" Guidewire

New 5F guide compatible Pronto LP 5F (Model 5015)

- Unique low-profile design offers 5F guide compatibility
- U-channel shaft for maximum aspiration



DISTAL SHAFT



PROXIMAL SHAFT

shown in 5F guide

0.014" Guidewire