



RUSCH® CRYSTALCLEAR™ PLUS TRACHEOSTOMY SET

Better comfort and increased safety for the patient and healthcare provider.

STANDARD TRACHEOSTOMY TUBE


Better comfort and increased safety for the patient and healthcare provider.

CRYSTALCLEAR™ PLUS TRACHEOSTOMY SET: ADVANCED TUBE SYSTEM WITH A FLEXIBLE INNER CANNULA MADE OF THERMOSENSITIVE PVC

The thermosensitive, soft plastic material adapts to the patient's anatomy. Secretions are removed simply and efficiently at an early stage of accumulation through frequent changes of the inner cannula, while protecting the inner tube lumen against contamination. Including an open and transparent fixation flange and an integrated rotating connector, this system is designed to provide better comfort and increased safety during use. The tubes feature low-pressure cuffs and come in intermediate sizes.


THE TUBE: FEATURES A LOW PRESSURE CUFF

- Thermosensitive tube material for greater flexibility
- Continuous X-ray marker along the tube shaft
- Integrated rotating 15 mm connector
- Soft connector. Open and transparent fixation flange for comfortable fit.
- Sterile
- Product is not made with natural rubber
- Single-use system to reduce the risk of cross-contamination

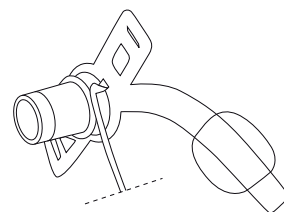


FLEXIBLE INNER CANNULA SYSTEM: COLORED STOPPER


- Flexible inner cannula system for ease of handling
- Color-coded stopper for quick and safe size identification
- Mucus/secretion can be collected inside the inner cannula, while protecting the inner tube lumen against contamination
- Supports infection prevention measures by preventing germs to grow into the tracheostomy tube by keeping the inner tube clear of secretions

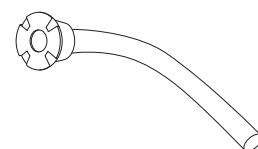

INSERTION AID: (OBTURATOR)

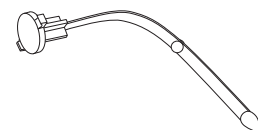

SOFT TRACHEOSTOMY TUBE NECKBAND:

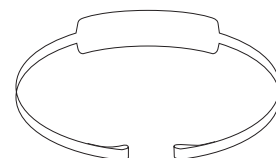

RUSCH® CRYSTALCLEAR™ PLUS TRACHEOSTOMY SET

Improves safety and supports infection prevention





SAFETY & EASE OF HANDLING

Our CrystalClear Plus tracheostomy tube set is a single-use system which helps to reduce the risk of infection, as the specially designed inner cannula system protects the inner tube lumen from secretion. Removal of secretions is simple and efficient at an early stage of accumulation through frequent changes of the inner cannula, and well-aligned accessory parts assist the patient in clearing mucus.

COMFORT & EFFICIENCY

Our CrystalClear Plus tracheostomy tube set is an advanced tube system with a flexible inner cannula system. It is distinguished by top quality, cost effectiveness and great wearing comfort. Simple and safe to use it renders the care of tracheostomy patients comfortable and efficient.

BROAD SIZE RANGE

Our CrystalClear Plus is available from sizes I.D. 6.0 to 10.0 mm, in 0.5 mm steps.



For more information visit
www.TELEFLEX.com



CRYSTALCLEAR PLUS HOSPITAL SET

with cuff



SET COMPONENTS

- 1 tracheostomy tube
- 1 flexible inner cannulas with colour-coded stopper
- soft neck band
- insertion aid (obturator)

HOSPITAL SET	ref.	ORDER size/i.D.	Qty
	121611	6.0–10.0	5
HOSPITAL SET INNER CANNULA	ref.	ORDER size/i.D.	Qty
	121310	6.0–10.0	10

CRYSTALCLEAR PLUS HOME CARE SET

with cuff



SET COMPONENTS

- 1 tracheostomy tube
- 3 flexible inner cannulas with colour-coded stopper
- soft neck band
- insertion aid (obturator)

HOME CARE SET	ref.	ORDER size/i.D.	Qty
	851611	6.0–10.0	1
HOME CARE SET INNER CANNULA	ref.	ORDER size/i.D.	Qty
	858610	6.0–10.0	2

ORDERING INFORMATION

CRYSTALCLEAR PLUS TRACHEOSTOMY TUBE

RUSCH

Adult	ORDER SIZE / I.D.	o.d.	total length (CLL=a+b+c)	DIM. A	DIM. B	DIM. C	Angle θ (approx.)	cuff ø *	inner cannula
	6.0	8.7 mm	64 mm	25.5 mm	0	38.5 mm	95°	19 mm	6.0
	6.5	8.7 mm	64 mm	25.5 mm	0	38.5 mm	95°	19 mm	6.5
	7.0	10.0 mm	70 mm	29.2 mm	0	40.8 mm	95°	21 mm	7.0
	7.5	10.7 mm	75.5 mm	30 mm	0	45.5 mm	95°	25 mm	7.5
	8.0	11.3 mm	81 mm	35.5 mm	0	45.5 mm	95°	28 mm	8.0
	8.5	11.3 mm	81 mm	35.5 mm	0	45.5 mm	95°	28 mm	8.5
	9.0	12.7 mm	87 mm	39.9 mm	0	47.1 mm	95°	28 mm	9.0
	9.5	12.7 mm	87 mm	39.9 mm	0	47.1 mm	95°	28 mm	9.5
	10.0	14.0 mm	95.5 mm	48.4 mm	0	47.1 mm	95°	30 mm	10.0

Teleflex, CrystalClear and Rusch are trademarks of Teleflex Incorporated or its affiliates, registered in the U.S. and/or other countries.

Teleflex is a global provider of medical products designed to enable healthcare providers to protect against infections and improve patient and provider safety. The company specializes in products and services for vascular access, respiratory, general and regional anesthesia, cardiac care, urology and surgery. Teleflex also provides specialty products for device manufacturers.

© 2015 Teleflex Incorporated. All rights reserved. LA-MC-000595



TELEFLEX LATIN AMERICA 3015 Carrington Mill Boulevard, Morrisville, NC 27560 USA
For product information please contact la.cs@teleflex.com or +1.919.433.4999
TELEFLEX.COM