

PRONTO® ***V3***
extraction catheter

DELIVERABILITY IN DESIGN



PRONTO V3

extraction catheter

Deliverability and Vessel Protection

The Pronto V3 was specially designed for deliverability during percutaneous thrombus aspiration with a self-centering tip to help prevent vessel wall adhesion. The full-length braided shaft enhances pushability and kink resistance while maintaining flexibility. This combination offers impressive deliverability when and where it's needed most.

Patented, Silva Tip

protects the vessel wall during introduction and aspiration

PTFE Lined Wire Lumen

for smooth guidewire movement

Shaft Construction

Proximal Shaft

Distal Shaft

Full-length Braided Shaft


Hydrophilic Coating

on the distal segment promotes smooth vessel navigation and trackability

Designed for Rapid Thrombus Removal

Uniformly Large Extraction Lumen

maximizes thrombus extraction force with no restrictions

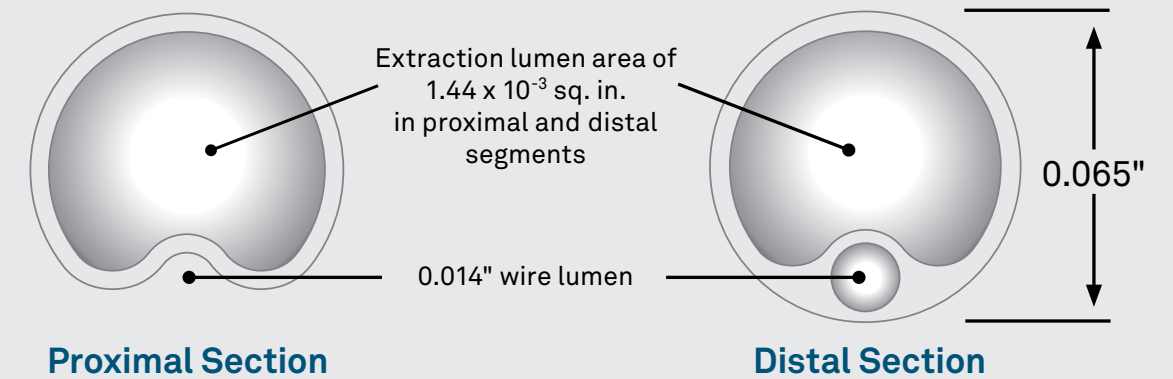


Aspiration Speed*

Pronto® V3	Export® Advance	PriorityOne™ AC	Export® AP	Xpress-Way™ RX
1.78cc/s	1.73cc/s	1.64cc/s	1.55cc/s	1.49cc/s

*Testing completed by Vascular Solutions, Inc. Data on file. Bench test results are not necessarily indicative of clinical performance.

Tapered, Silva Tip

is designed to improve navigation and self-centers for vessel protection during delivery

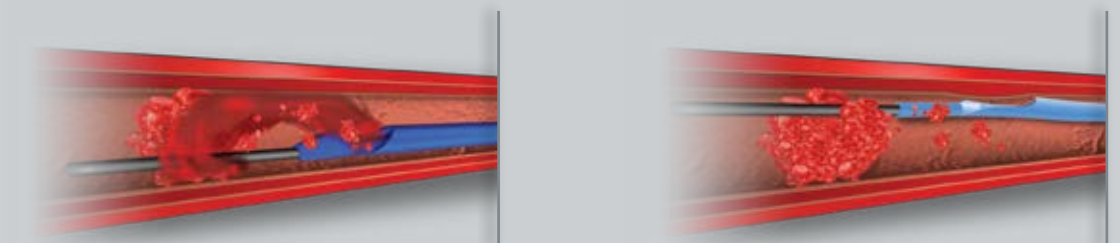

Pronto Silva tip

Standard open-end extraction catheter tip

PRONTO[®] V3

extraction catheter

The Pronto V3 extraction catheter is indicated for the removal of fresh, soft emboli and thrombi from vessels in the coronary and peripheral vasculature.

Model	Description	Guide Catheter Compatibility	Guidewire Compatibility	Working Length	Rapid Exchange Section	Box Quantity
5003	Pronto V3 Extraction Catheter	≥6F (0.070" min. I.D.)	≤0.014"	140cm	20cm	1

Packaged in quantities of 1 unit per box.

Please see the Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions.

CAUTION: Federal Law (U.S.A.) restricts this device to sale by or on the order of a physician.

Pronto is a registered trademark of Vascular Solutions, Inc. All other trademarks and registered trademarks are property of their respective owners.

©2016 Vascular Solutions, Inc. All rights reserved. ML1559 Rev. K 05/16



Vascular Solutions, Inc.
6464 Sycamore Court North
Minneapolis, Minnesota 55369 USA

Customer Service:
United States: 888.240.6001
International: (001) 763.656.4298
customerservice@vasc.com

www.vasc.com