

Ordering Information

Arrow Simplicity® Color-Coded Micro-Introducer

PRODUCT NUMBER	DESCRIPTION	QUANTITY	SETS/CASE
MP-04450-AS	Color Coded Micro-Introducer Set	5 ea. per box	1

Arrow Edge® Hemodialysis Catheterization Set and Arrow Simplicity Combination Product

COLOR CODE	PRODUCT NUMBER	DESCRIPTION	TIP-TO-CUFF INSERTION LENGTH	SETS/CASE
19	CS-15242-IM	Arrow 15 Fr. x 9.5" (24 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	19 cm	1
23	CS-15282-IM	Arrow 15 Fr. x 11" (28 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	23 cm	1
27	CS-15322-IM	Arrow 15 Fr. x 12.5" (32 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	27 cm	1
31	CS-15362-IM	Arrow 15 Fr. x 14" (36 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	31 cm	1

Arrow® Cannon® II Plus Hemodialysis Catheterization Set and Arrow Simplicity Combination Product

COLOR CODE	PRODUCT NUMBER	DESCRIPTION	TIP-TO-CUFF INSERTION LENGTH	SETS/CASE
19	CS-15242-SPM	Arrow 15 Fr. x 9.5" (24 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	19 cm	1
23	CS-15282-SPM	Arrow 15 Fr. x 11" (28 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	23 cm	1
27	CS-15322-SPM	Arrow 15 Fr. x 12.5" (32 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	27 cm	1
31	CS-15362-SPM	Arrow 15 Fr. x 14" (36 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	31 cm	1

Arrow® NextStep® Antegrade Hemodialysis Catheterization Set and Arrow Simplicity Combination Product

COLOR CODE	PRODUCT NUMBER	DESCRIPTION	TIP-TO-CUFF INSERTION LENGTH	SETS/CASE
19	CS-15192-SFXM	Arrow 15 Fr. x 7.5" (19 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	19 cm	1
23	CS-15232-SFXM	Arrow 15 Fr. x 9" (23 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	23 cm	1
27	CS-15272-SFXM	Arrow 15 Fr. x 10.5" (27 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	27 cm	1
31	CS-15312-SFXM	Arrow 15 Fr. x 12" (31 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	31 cm	1

Arrow® NextStep® Retrograde Hemodialysis Catheterization Set and Arrow Simplicity Combination Product

COLOR CODE	PRODUCT NUMBER	DESCRIPTION	TIP-TO-CUFF INSERTION LENGTH	SETS/CASE
19	CS-15192-XM	Arrow 15 Fr. x 7.5" (19 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	19 cm	1
23	CS-15232-XM	Arrow 15 Fr. x 9" (23 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	23 cm	1
27	CS-15272-XM	Arrow 15 Fr. x 10.5" (27 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	27 cm	1
31	CS-15312-XM	Arrow 15 Fr. x 12" (31 cm) radiopaque polyurethane catheter 21 Ga. Micro-Introducer set	31 cm	1

Caution: Federal law restricts this device to sale by or on the order of a physician. Contents of unopened, undamaged package are sterile. Disposable. Refer to package insert for current warnings, indications, contraindications, precautions and Instructions For Use.

Ordering Information (cont.)

Arrow-Clark™ VectorFlow® Chronic Hemodialysis Catheter and Arrow Simplicity

COLOR CODE	PRODUCT NUMBER	DESCRIPTION	TIP-TO-CUFF INSERTION LENGTH	SETS/CASE
19	CS-15192-VFIM	VF -Std w/ Simplicity HEMODIALYSIS KIT: 2-L 15 FR X 19 CM ANTEGRADE with Simplicity	19 cm	1
23	CS-15232-VFIM	VF -Std w/ Simplicity HEMODIALYSIS KIT: 2-L 15 FR X 23 CM ANTEGRADE with Simplicity	23 cm	1
27	CS-15272-VFIM	VF -Std w/ Simplicity HEMODIALYSIS KIT: 2-L 15 FR X 27 CM ANTEGRADE with Simplicity	27 cm	1
31	CS-15312-VFIM	VF -Std w/ Simplicity HEMODIALYSIS KIT: 2-L 15 FR X 31 CM ANTEGRADE with Simplicity	31 cm	1
19	CS-15192-VFM	HEMODIALYSIS KIT: 2-L 15 FR X 19 CM RETROADE with Simplicity	19 cm	1
23	CS-15232-VFM	HEMODIALYSIS KIT: 2-L 15 FR X 23 CM RETROADE with Simplicity	23 cm	1
27	CS-15272-VFM	HEMODIALYSIS KIT: 2-L 15 FR X 27 CM RETROADE with Simplicity	27 cm	1
31	CS-15312-VFM	HEMODIALYSIS KIT: 2-L 15 FR X 31 CM RETROADE with Simplicity	31 cm	1

References:

1. Refer to Arrow chronic hemodialysis catheter package insert*

Teleflex is a global provider of medical technologies designed to improve the health and quality of people's lives. We apply purpose-driven innovation – a relentless pursuit of identifying unmet clinical needs – to benefit patients and healthcare providers. Our portfolio is diverse, with solutions in the fields of vascular and interventional access, surgical, anesthesia, cardiac care, urology, emergency medicine and respiratory care. Teleflex employees worldwide are united in the understanding that what we do every day makes a difference. For more information, please visit teleflex.com.

Teleflex is the home of Arrow®, Deknatel®, Hudson RCI®, LMA®, Pilling®, Pronto, Rüschi®, UroLift®, and Weck® – trusted brands united by a common sense of purpose.


Please see the Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions

CAUTION: Federal (USA) law restricts this device to sale by or on the order of a physician.]

Teleflex, the Teleflex logo, Arrow, Arrow Simplicity, Arrow Edge, Cannon and NextStep are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries.

Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative. Revised: 12/2019

© 2019 Teleflex Incorporated. All rights reserved. MC-005149 * Rev 1

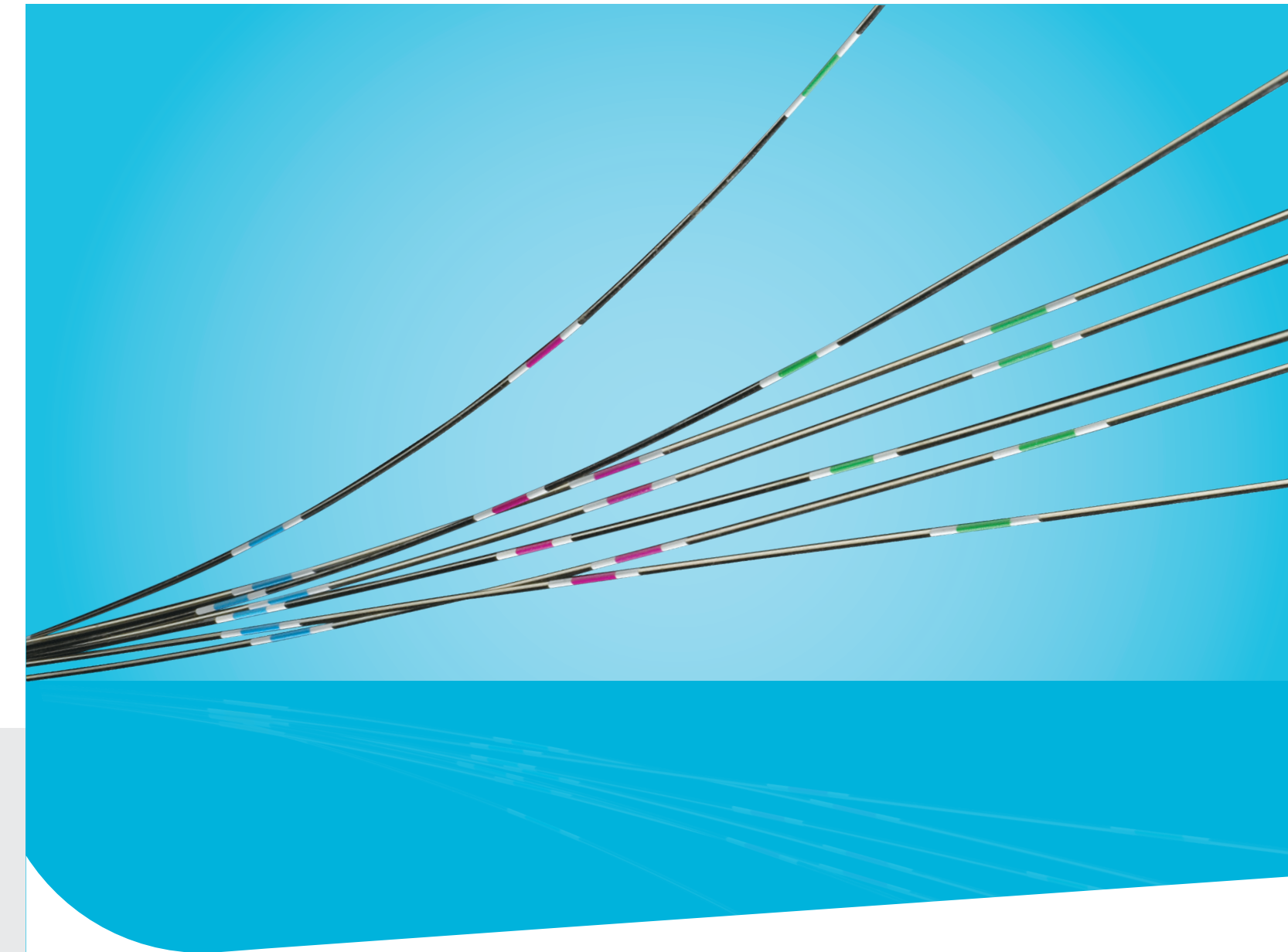

Arrow Simplicity® Color-Coded Micro-Introducer Set Catheterization made simpler




Arrow Simplicity® Color-Coded Micro-Introducer

Designed to speed catheterization

The Arrow Simplicity Color-Coded Micro-Introducer Set streamlines catheterization by facilitating optimal catheter positioning and standardizing catheter length selection.

Enhances visibility during ultrasound guided insertions
21 Ga. x 2.75" thin wall echogenic needle with green hub and bevel indicator.

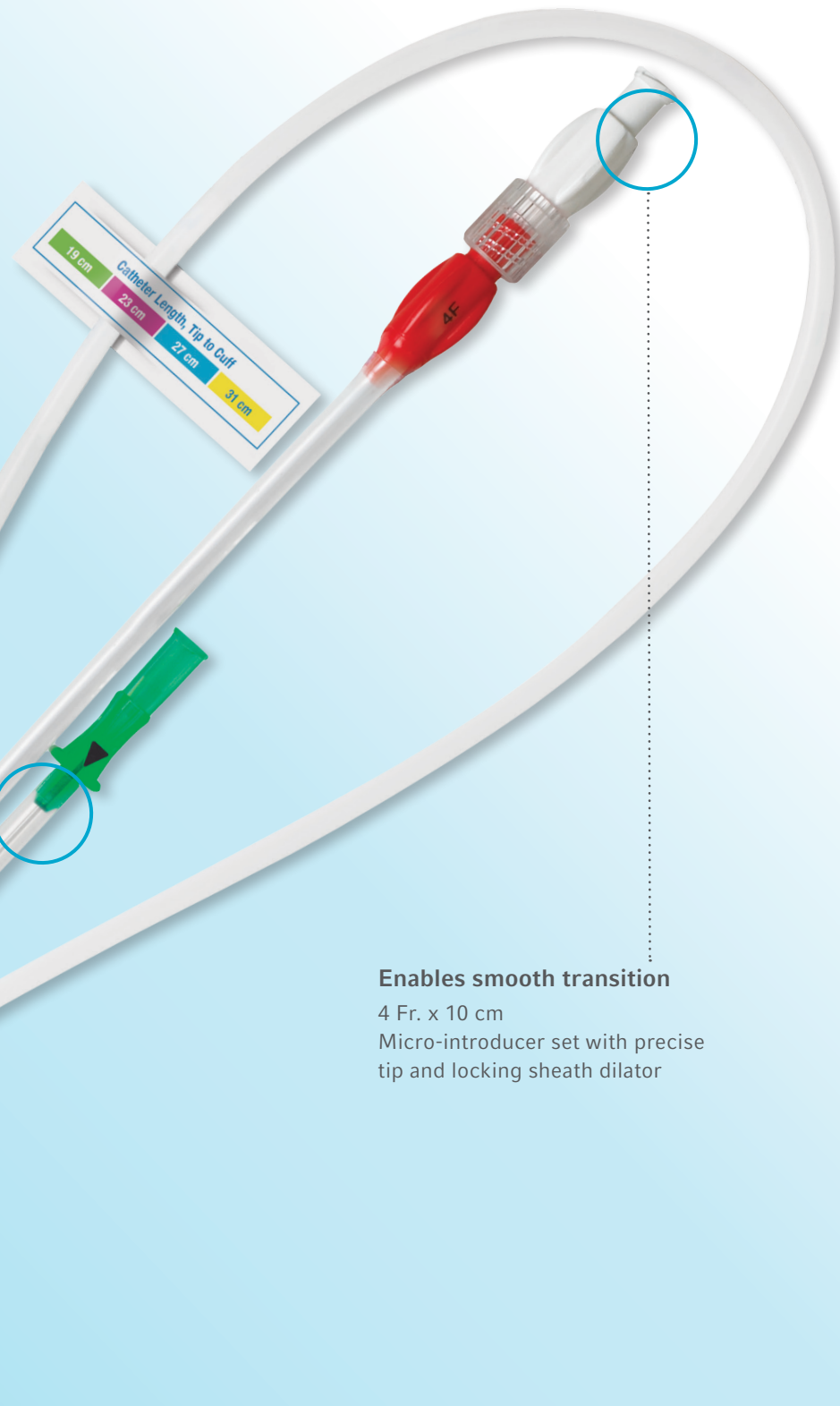

Enables smooth transition
4 Fr. x 10 cm
Micro-introducer set with precise tip and locking sheath dilator

Arrow Simplicity® Color-Coded Micro-Introducer



Positioning made simpler

With Arrow Simplicity Color-Coded Micro-Introducer you can determine the optimal catheter length required in three easy steps.

STEP 1: Insert floppy coiled end of Arrow Simplicity color-coded micro-introducer guidewire distal tip to desired location.¹


STEP 2: Position Arrow Simplicity color-coded micro-introducer guidewire over patient's chest to desired tunnel tract configuration.

STEP 3: Determine the length of Arrow chronic hemodialysis catheter you require by matching the color band on the Arrow Simplicity guidewire at the exit site with the color-coded Arrow chronic hemodialysis catheter package.

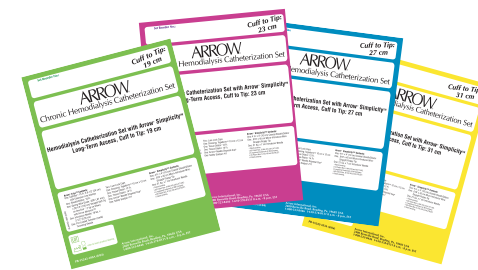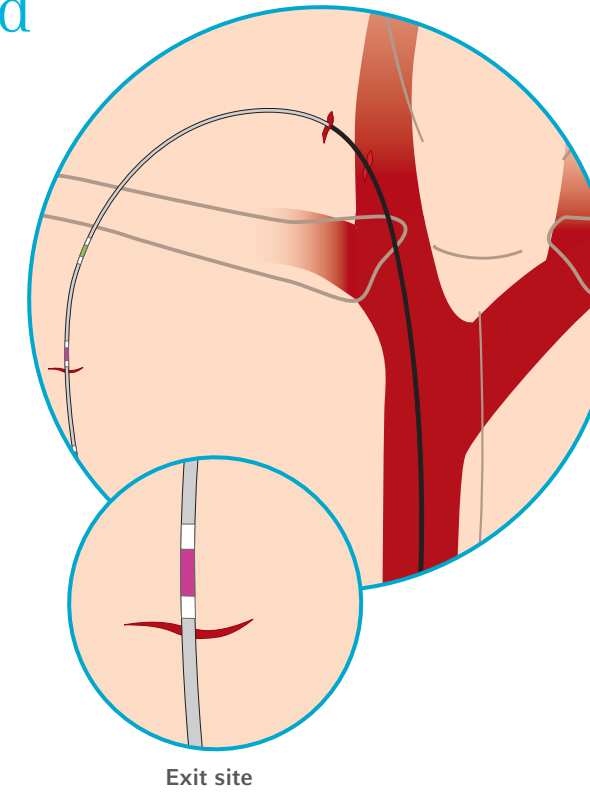

Colored bands on guidewire correspond to color-coded catheter packaging
Arrow Simplicity can only be used with Arrow chronic hemodialysis catheters

Heightens visibility
0.018" x 45 cm stainless steel guidewire with soft, flexible, tapered platinum coil tip

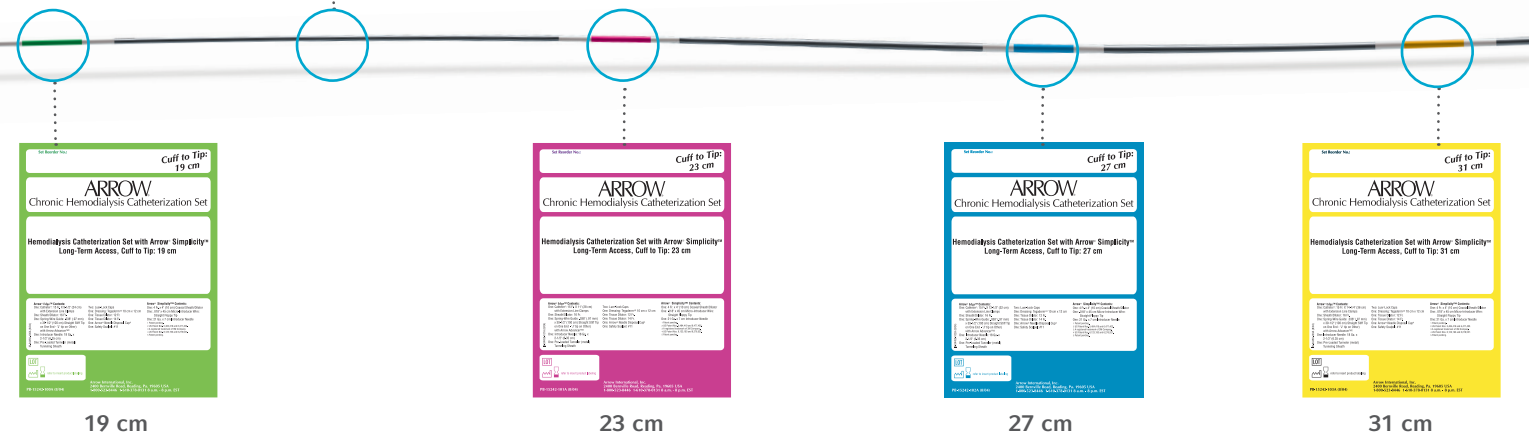

Arrow Simplicity® Micro-Introducer Set

Contents:

- 1 Introducer Needle: 21 Ga. x 2.75" (7 cm) thin wall echo tip needle with bevel indicator and green hub
- 1 Coaxial Sheath/Dilator: 4 Fr. x 4" (10 cm)
- 1 Micro-Introducer Wire Guide: 0.018" (0.46 mm) x 45 cm stainless steel guidewire with platinum coil and color band markings



Arrow Simplicity Micro-Introducer Set
MP-04450-AS

To learn more about the Arrow® Simplicity Color-Coded Micro-Introducer Set, contact your Teleflex representative.