

Rusch®
EquipLite™
 Standard Disposable
 Laryngoscope Blade

Minimize Risk, Maximize Outcomes



Rusch® EquipLite™ Blades are ideal for O.R. intubations, crash carts and emergency situations

Call 1.866.246.6990 to place an order

A Bright Solution to Reduce the Risk of Cross-Contamination

Intubation of the trachea has been a risky cross-contamination procedure over the past decade because no perfect decontamination procedures exist.¹ Single-use devices offer the opportunity to reduce potential for infection associated with reusable devices.

The Rusch® EquipLite™ Disposable Metal Laryngoscope Blade from Teleflex offers users the performance and feel of a reusable metal blade.

The Rusch® EquipLite™ Blade eliminates the hassles associated with reusable blades, including the time and cost of reprocessing, without sacrificing strength or light quality.

The Performance of Metal with the Convenience of a Single-Use Blade

Built with solid metal construction, the Rusch® EquipLite™ Disposable Standard Blade offers outstanding performance. The blade provides a simple, reliable click locking and unlocking mechanism and features an integrated LED bulb that provides optimal light output.

Features for Confident and Convenient Use

The Rusch® EquipLite™ Laryngoscope Blade is constructed of surgical grade stainless steel with a metal reinforced heel and a beveled blade tip. The blade comes in a large range of color-coded sizes and is compatible with standard/conventional handles.

The single-use blades have a medium profile and height and feature LED lighting technology that provides strong, focused and reliable illumination while reducing battery consumption.



BENEFITS



The Clinician
 Design helps inspire clinician confidence during critical procedures




Your Institution
 Cost-neutral* way to combat intubation infection risk



The Patient
 Promotes positive patient outcomes

**Individual costs and savings may vary based on a number of factors, including specific markets, services and rates.*

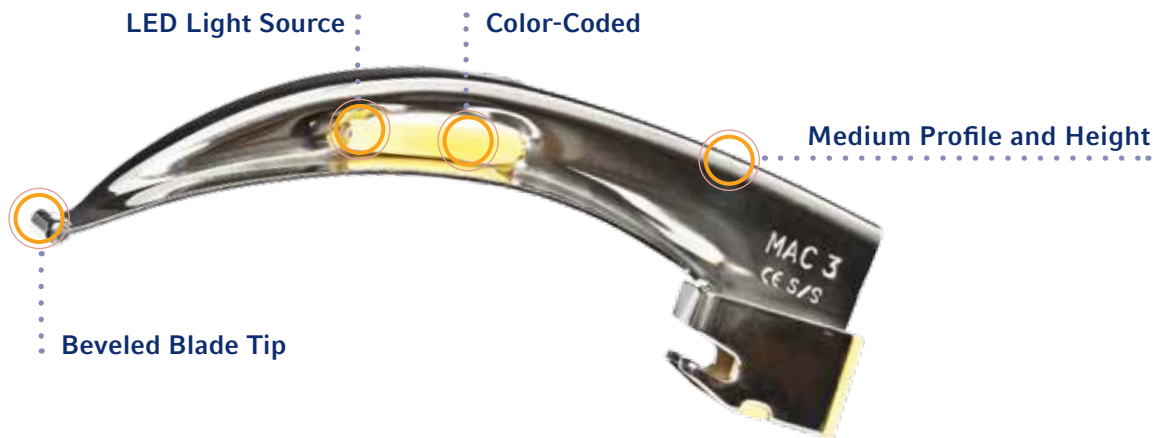

Rusch® EquipLite™ Fast Facts

Compatible with standard/conventional handles. Conforms to ISO standards to allow use with compatible handles

The bright blue-white LED light source enhances tissue differentiation and since LED is a cool light source, it eliminates the risk of burning a patient's mouth. The LED does not contain a bulb filament, allowing for a uniform light with no shadowing

Metal reinforced heel "clicks" into place, ensuring the blade is properly seated and provides additional strength at this stress point

Medium profile and height allows better maneuverability and visualization, decreases the chances of knocking teeth, and improves handling when facing limited mouth openings

Beveled blade tip reduces potential trauma to the patient

Macintosh Blades

MEDICAL P/N	DESCRIPTION	OVERALL LENGTH (MM)	DISTAL WIDTH (MM)	CASE QTY	COLOR-CODE
004651000	Rusch® EquipLite™ Mac 0	77	9.0	20	Peach
004651001	Rusch® EquipLite™ Mac 1	92	9.0	20	White
004651002	Rusch® EquipLite™ Mac 2	101	12.0	20	Turquoise
004651003	Rusch® EquipLite™ Mac 3	131	15.0	20	Yellow
004651035	Rusch® EquipLite™ Mac 3.5	147	18.5	20	Coral
004651004	Rusch® EquipLite™ Mac 4	155	15.0	20	Pink
004651005	Rusch® EquipLite™ Mac 5	172	15.0	20	Dark Yellow

Miller Blades

MEDICAL P/N	DESCRIPTION	OVERALL LENGTH (MM)	DISTAL WIDTH (MM)	CASE QTY	COLOR-CODE
004650010	Rusch® EquipLite™ Miller 00	64	11.5	20	Beige
004650000	Rusch® EquipLite™ Miller 0	75	11.5	20	Purple
004650001	Rusch® EquipLite™ Miller 1	102	11.5	20	Orange
004650002	Rusch® EquipLite™ Miller 2	155	14.0	20	Gray
004650003	Rusch® EquipLite™ Miller 3	195	14.0	20	Green
004650004	Rusch® EquipLite™ Miller 4	205	14.0	20	Lavender

REFERENCE

1 Disposable laryngoscope blades do not interfere with ease of intubation in scheduled general anesthesia patients. *European Journal of Anaesthesiology* 2003;20: 731-735.



TELEFLEX.COM/LARYNGOSCOPES

EquipLite is a trademark of Truphatek International Ltd. Teleflex and Rusch are trademarks or registered trademarks of Teleflex Incorporated or its affiliates. © 2015 Teleflex Incorporated. All rights reserved. MC-000286