



532221
Pilling™ Adult Laser
Dedo Laryngoscope

Pilling™
Vitality Dedo Laser Laryngoscopes
 The Next Generation of Pilling™ Laryngoscopes are Here

The Pilling™ Vitality line of laryngoscopes offers your familiar Dedo access with the same proximal opening, length, and internal diameter dimensions while providing an improved accessory interface and main tube construction.

Features and Benefits

- Modular accessories
- Color Coding
- All scopes are laser compatible
- All scopes have jet ventilation port
- Solid stainless steel construction




532233
Pilling™ Ossoff-Karlan Posterior
Commissure Laryngoscope



532224
Pilling™ Adolescent
Laser Dedo Laryngoscope



Accessories are flush to the handle to provide an ergonomic grip and keep cables and supply tubing out of the way

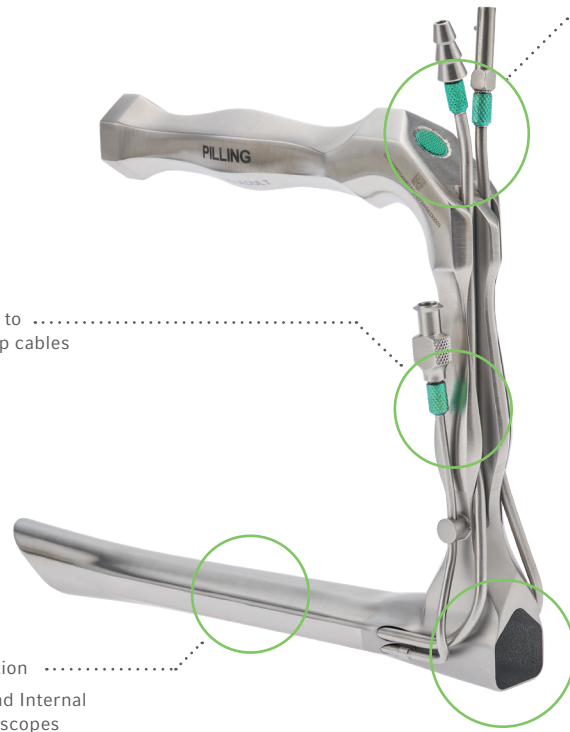

Color coded accessories for better organization and scope assembly



Stainless steel main tube construction
Main tube provides same access and Internal diameter as legacy Pilling™ laryngoscopes

Three accessory channels integrated into the main tube for light, smoke evacuation, and jet ventilation
Clear visual field with all accessories inserted

NEW Pilling™ Vitality Dedo and Ossoff-Karlan Laryngoscopes

| ITEM NO. | DESCRIPTION | PREVIOUS ITEM NO. | OVERALL LENGTH | MIN I.D. | NOMINAL PROXIMAL DIMENSIONS | MAX O.D. |
|----------|---|-------------------------|----------------|----------|-----------------------------|----------|
| 532221 | Adult Dedo Laser Micro Laryngoscope | 522221, 522225P, 522228 | 178 mm | 10.7 mm | 19.0 mm x 19.0 mm | 28.4 mm |
| 532222 | Child Dedo Micro Laser Laryngoscope | 522222 | 119 mm | 7.1 mm | 15.0 mm x 16.0mm | 23.5 mm |
| 532224 | Adolescent Dedo Micro Laser Laryngoscope | 522224 | 145 mm | 9.1 mm | 18.0 mm x 17.0 mm | 24.6 mm |
| 532237 | Child Ossoff-Karlan Laser Laryngoscope | 522237 | 130 mm | 9.8 mm | 18.0 mm x 18.0 mm | 24.3 mm |
| 532238 | Adolescent Ossoff-Karlan Laser Laryngoscope | 522238 | 152 mm | 9.7 mm | 21.0 mm x 20.0 mm | 27.1 mm |
| 532231 | Adult Ossoff-Karlan Laser Laryngoscope | 522231 | 178 mm | 11.9 mm | 22.0 mm x 20.0 mm | 30.9 mm |
| 532233 | Ossoff-Karlan Posterior Commissure Laryngoscope | 522233 | 177 mm | 12.6 mm | 22.0 mm x 19.0 mm | 31.8 mm |
| 532243 | Dedo Laser Extra-Long Laryngoscope | 522243 | 197 mm | 12.6 mm | 20.0 mm x 20.0 mm | 28.6 mm |

Individual Components

| SCOPE ONLY | LIGHT CARRIER | SMOKE EVACUATION TUBE | JET VENTILATION | COLOR CODE |
|------------|---------------|-----------------------|-----------------|------------|
| 532221S | 537525 | 538525 | 539525 | ● Green |
| 532222S | 537524 | 538524 | 539524 | ● Pink |
| 532224S | | | | |
| 532237S | 537524 | 538524 | 539524 | ● Pink |
| 532238S | | | | |
| 532231S | 537525 | 538525 | 539525 | ● Green |
| 532233S | 537525 | | | |
| 532243S | 537525 | | | |



To view NEW Pilling™ Vitality Laryngoscope products go to: myTeleflex.com

Teleflex, the Teleflex logo, and Pilling are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Information in this document is not a substitute for the product Instructions for Use. The products in this document may not be available in all countries. Please contact your local representative. All data current at time of printing (01/2024). Subject to technical changes without further notice.
© 2024 Teleflex Incorporated. All rights reserved. MC-009490