



ARROWGARD⁺ BLUE[®] MAC[™]

In the ErgoPack[™] Maximal Barrier Precautions Procedural Tray

RISK REDUCTION

Infection Protection Features:

- ARROWg⁺ard Blue MAC
- Procedural Stop Sign
- Central Line Insertion Checklist
- Central Line Insertion Practices Adherence Monitoring Form (CLIP Sheet)
- Hi-Lite Orange[®] ChloraPrep[®]
- Alcohol Gel Hand Sanitizer
- Sterile Hand Towel
- Maximal Barrier Set
 - Surgical Cap
 - Mask with Face Shield
 - Surgical Gown
 - Exclusive Easy-to-use Head-to-toe Drape
- BioPatch[®]
- Tegaderm[®] Semi-permeable Transparent Dressing
- HemoHopper[®] Fluid Receptacle
- Cath-Gard[®] Contamination Shield

Sharps Safety Features:

- SharpsAway[™] II Locking Disposal Cup
- SafetyGlide[®] Injection Needles
- Safety Scalpel

The ARROWg⁺ard Blue MAC ErgoPack provides 30 essential tools systematically organized for procedural ease. The ErgoPack integrates the latest infection prevention and safety measures recommended by the CDC, SHEA/IDSA, IHI, JCHAO, NHSN and OSHA into one convenient bundle.

The ARROWg⁺ard Blue MAC is a unique product designed to meet the needs of a variety of patient types — from trauma to cardiac. The catheter combines the access of a sheath introducer with the high-flow

lumens of a central line. The MAC's hemostasis valve allows for easy access for additional devices, such as a Thermodilution catheter or MAC Companion catheter adding up to 3 additional lumens. The 9 French distal lumen accommodates large fluid volumes when needed, while the 12 gauge proximal lumen provides a separate access for incompatibles fluid administration.

The ARROWg⁺ard Blue MAC catheter uses Teleflex's proprietary chlorhexidine based antimicrobial technology proven to protect your patients from the known pathogens of Catheter-Related Bloodstream Infections.




ARROWg⁺ard Blue Multi-Lumen Access Catheter (MAC)



ErgoPack Maximal Barrier Precautions Procedural Tray

ORDERING INFORMATION


ERGOPACK MAC PRODUCTS*

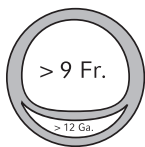

Product Number	Description	Per Case
CDC-21242-X1A	9 Fr. x 10 cm. 2-Lumen with Integral Hemostasis Valve	5
<i>For Sale Outside the US:</i>		
CDA-21242-1A	9 Fr. x 10 cm. 2-Lumen with Integral Hemostasis Valve	5

*For a customized option through our Arrow Select Kits program, contact your local Teleflex Vascular Access Sales Specialist.

MAC COMPANION PRODUCTS AND ACCESSORIES

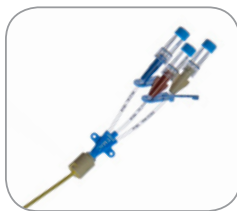

Product Number	Description	Fr.	Length	Per Case
ML-00701	Single-Lumen; 12 ga.	7 Fr.	20 cm.	5
ML-00702	Double-Lumen; 16, 16 ga.	7 Fr.	20 cm.	5
ML-00703	Triple-Lumen; 16, 18, 18 ga.	7 Fr.	20 cm.	5
AO-09000	Hemostasis Valve Oburator Cap			50

CATH-GARD CONTAMINATION SHIELD



Product Number	For use with Catheter:	Length	Connection	Color	Per Case
ST-09807	7-7.5 Fr.	80 cm.	Tuohy-Borst	Blue	25
ST-09880	7-7.5 Fr.	80 cm.	No lock	Blue	25
ST-09870	7 Fr.	80 cm.	TwistLock	Blue	25
ST-09875	7.5-8 Fr.	80 cm.	TwistLock	Yellow	25
ST-09885	8-8.5 Fr.	80 cm.	TwistLock	Black	25
ST-09847	4-7 Fr.	80 cm.	Tuohy-Borst	Blue	10

TO PLACE AN ORDER CONTACT CUSTOMER SERVICE AT 1.866.246.6990

Caution: US federal law restricts these devices to sale by or on the order of a physician.

TELEFLEX MEDICAL Research Triangle Park, NC 27709
Phone: 866.246.6990 Fax: 800.399.1028 Intl: 919.361.3961

TELEFLEXMEDICAL.COM

BioPatch is a registered trademark of Johnson and Johnson Corporation. Hi-Lite Orange and ChloroPrep are registered trademarks of CareFusion. SafetyGlide is a registered trademark of Becton, Dickinson, and Company. Tegaderm is a registered trademark of 3M Company. Arrow, ARROWgard Blue, Cath-Gard, ErgoPack, HemoHopper, MAC, SharpsAway and Teleflex are trademarks or registered trademarks of Teleflex Incorporated or its affiliates.

©2010 Teleflex Incorporated. All rights reserved. 2010-0376 v2

Teleflex®